

Multilingual manuals



Scan QR code for supported language manuals (Google recommended)

QR-Code scannen für Anleitungen in unterstützten Sprachen (Google empfohlen)

Escanea el código QR para manuales en idiomas admitidos (se recomienda Google)

Zeskanuj kod QR, aby uzyskać instrukcje w obsługiwanych językach (zalecane Google)

Desteklenen dillerdeki kılavuzlar için QR kodu tara (Google önerilir)

Скануйте QR-код для інструкцій підтримуваними мовами (рекомендовано Google)

Imbas kod QR untuk manual dalam bahasa yang disokong (disyorkan Google)

Pindai QR code untuk panduan dalam bahasa yang didukung (direkomendasikan Google)

I-scan ang QR code para sa mga manual sa tinutukoy na wika (inirerekomenda ang Google)

Quét mã QR để xem hướng dẫn bằng ngôn ngữ được hỗ trợ (khuyến dùng Google)

สแกนรหัส QR เพื่อรับคู่มือภาษาที่รองรับ (แนะนำ Google)

扫码获取支持语言说明书 (推荐Google)

Enjoy the Power-On Experience



Foot cover zipper
(Removable for cleaning)

Adapter interface
Connect the adapter here

Message mode switch button
Touch this button to switch
message modes
(three modes in total).

Hot compress button
Touch this button to
turn on / off heat

Airbag gear switch button
Touch this button to switch
the airbag gear
(three gear in total)

Timer button
Touch this button to set the massage time
(15 min/30 min)

Controller

Foot Covers
Place your feet inside the
foot covers and
begin your massage experience

* The schematic diagrams are for reference only. Actual products may vary. The actual products shall prevail.

How to use

Foot massage 1–2 times a day for 15–30 minutes each time can relieve muscle pain, promote blood circulation, and improve insomnia. It is recommended not to use it for more than 30 minutes at a time, as prolonged use may cause slight pain.

- After connecting the power supply, the power indicator light will flash and enter standby mode.
- Remove the touch panel protective film and press the power button to start. The default air pressure intensity is low (L), the default mode is mode 1, and it automatically turns off after 15 minutes (red light).
- Press the "Timer" button to switch to 30-minute mode.
- Find the most suitable function combination through the "Massage Mode" and "Air Pressure Intensity" buttons (L/M/H/OFF).
- Press the "Heat" button to turn on or off the heating function.
Note: The heating function must be used together with the massage function and cannot be used alone.
Tip: Pressing the button with your thumb can increase touch sensitivity.
- The device is compatible with men's shoes up to US size 13 (EU size 47).
- The cloth cover is removable and washable (see picture for details).

Important Safety Instructions

Please read all instructions before use.

- This appliance is not intended for use by persons (including children) with limited physical, sensory or mental abilities, unless they are supervised or instructed by a person responsible for their safety.
- Children must be supervised at all times and are prohibited from using the device unattended.
- WARNING** — To reduce the risk of burns; fire; electric shock; or injury to persons:
a) Do not leave the device unattended when it is plugged in. Unplug the device when not in use and before installing or removing parts.
b) Do not operate the device under a blanket or pillow. Overheating may cause fire, electric shock or personal injury.
c) Children, patients or disabled people must be closely supervised when using or touching the device.
d) Use this device only as described in this manual. Do not use accessories not recommended by the manufacturer.
e) Stop using the device immediately if the power cord or plug is damaged, the device is abnormal, has been dropped or soaked in water.
f) Do not lift the device by the power cord.
g) Keep the power cord away from hot surfaces.
h) Make sure the device vents are not blocked and keep them free of lint, hair and other debris.
i) Do not throw foreign objects into the device openings.
j) Do not use outdoors.
k) Do not use in an environment where aerosol sprays or oxygen is administered.
l) To turn off the device, turn it to "off" before unplugging it.
- Technical Specification:** The device operates at safety extra-low voltage (see label for nominal voltage).
- Heated surfaces may pose a risk to people with heat insensitivity and should be used with caution.
- WARNING:** Hot surfaces may cause severe burns. Do not use on skin areas with reduced sensitivity or poor circulation. Children or incapacitated persons must be supervised when using the heating function.
- WARNING:** Do not stand on the device, use it in a sitting position only.

Please keep this manual properly

SKG Health Technologies Co., Ltd.

Address: 23A Floor, Building 3, Zhongke R&D Park, No. 009,
Gaoxin South 1st Road, High-tech Zone Community, Yuehai Street,
Nanshan District, Shenzhen, Guangdong, China.

Website : <https://global.skg.com/>



Made in China

SKG Foot Massager – Product Information

Specification

Model: YS100	Input: 12V ≈ 2.5A
Rated Power: 30W	Heating Area: Bottom
Time: 15 minutes / 30 minutes	Material: EVA + Leather
Heating Temperature: ≤55°C* (131°F)	Color: Dark blue

Package Contents

SKG YS100 Foot Massager ×1
Adapter ×1
Manual ×1
Remote Control*1

Function

Time: Choose 15 or 30 minutes with automatic shutdown.
Massage Mode: Low (L) / Medium (M) / High (H) / Off (OFF).
Air Pressure Intensity: Low (L) / Medium (M) / High (H) / Off (OFF).
Heat: Press to turn on / off the heat function.

Maintenance and care

- It is forbidden to use non-nominal voltage (DC 12 V / 2.5 A).
- Non-professionals are prohibited from repairing the device or replacing parts.
- If the power adapter is damaged, replace it only through the manufacturer's certified service.
- Do not use/store the device in a humid or high temperature environment.
- Please unplug the power plug after use.
- Limit continuous use to 30 minutes per session.
- Please stop using it immediately if you feel uncomfortable.
- Please turn off the power immediately when the device is hot or emits odor.
- Pregnant women, patients with heart disease, diabetes, or thrombosis, and individuals with pacemakers or metal implants should consult a doctor before use.
- Do not touch the device with sharp objects.
- Use a soft cloth and neutral detergent when cleaning.
- It is forbidden to stand or place heavy objects on the device.

Troubleshooting

Frequently Asked Questions	Reason	Solution
Device not working	The power supply is not connected / the plug is in poor contact	Re-plug the power cord and check the power switch.
Abnormal massage intensity	Improper foot position / incorrect intensity selection	Adjust your posture and start from low intensity
Heating function fails or is too slow	International safety standards limit heating speed	This is normal for safety. Allow extra heating time in cold environments
Touch panel malfunction	The protective film is not torn off. The fingers are wet / the touch area is not enough	Remove the protective film and keep your fingers dry

Product Disposal Precautions

Please support the recycling of rechargeable batteries (lithium-ion batteries)

	Danger	The battery used in this product is a dedicated rechargeable battery and cannot be used in other products. Please do not charge the rechargeable battery that has been removed. - Do not throw it into fire to overheat it - Do not pierce it with nails, impact it with gravity, or disassemble it - Do not connect the positive and negative poles of the battery directly with metal - Do not store it with necklaces, hairpins, etc. - Do not charge, use or place it near a stove, under the scorching sun, or in other high-temperature places. The above situations may cause heating, fire, rupture, etc.
	Warning	Please place the removed rechargeable battery out of reach of infants and young children to prevent accidental swallowing. - If swallowed by mistake, please consult a doctor immediately. - If the battery fluid leaks, do not touch the liquid with bare hands, and handle it as follows. - When the battery fluid enters the eyes, there is a risk of blindness. Do not rub the eyes directly with your hands. Immediately wash with clean water and consult a doctor. When the battery fluid gets on the body or clothes, it may cause skin inflammation or injury. Please rinse with clean water. If you still feel uncomfortable, consult a doctor immediately. - The product can only be disassembled when it is discarded. Failure to comply may cause fire, electric shock, injury and other accidents. If repair is required (except for battery replacement), please consult the store where you originally purchased this product .

The lithium-ion battery used in this product is a valuable resource that can be recycled. When you are finished using it, please hand it over to a professional organization that can recycle and recycle rechargeable batteries. Before discarding the device, the battery must be removed from the device. When removing the battery, the device must be powered off. The battery should be disposed of safely.

Removal steps:

- Unplug the charging cable and remove the product cover
- Remove the upper case and disconnect the battery power cord
- Remove the battery

CE This device was tested for typical body-worn operations. To comply with RF exposure requirements, a minimum separation distance of 0cm must be maintained between the user's body and the handset, including the antenna. Third-party belt-clips, holsters, and similar accessories used by this device should not contain any metallic components. Body-worn accessories that do not meet these requirements may not comply with RF exposure requirements and should be avoided. Use only the supplied or an approved antenna. This device in complies with the essential requirements and other relevant provisions of Directive 2014/53/EU. All essential radio test suites have been carried out.


1. CAUTION : RISK OF EXPLOSION IF BATTERY IS REPLACED BY AN INCORRECT TYPE. DISPOSE OF USED BATTERIES ACCORDING TO LOCAL REGULATIONS.

2. The device complies with RF specifications when the device used at 0mm form your body EU Compliance Statement:

SKG Health Technology Co., Ltd. Here by declares that this device is in compliance with essential requirements and other relevant provisions of the Directive 2014/53/EU. A copy of the EU Declaration of Conformity is available online at <https://www.skg.com/pages/skg-declaration-of-conformity>.

Old electrical appliances must not be disposed of together with the residual waste, but have to be disposed of separately. The disposal at the communal collecting point via private persons is for free.

The owner of old appliances is responsible to bring the appliances to these collecting points or to similar collection points. With this little personal effort, you contribute to recycle valuable raw materials and the treatment of toxic substances.

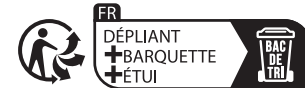

EU REP **SUCCESS COURIER SL**
CALLE RIO TORMES NUM.1, PLANTA 1, DERECHA,
OFICINA 3, Fuenlabrada, Madrid, 28947 Spain
successservice2@hotmail.com

OST_EU_2023021400009