




HealthyWhite

PHILIPS
sonicare

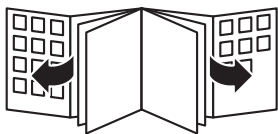1
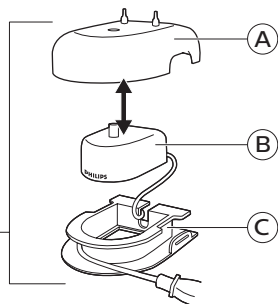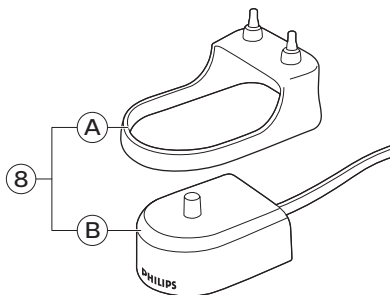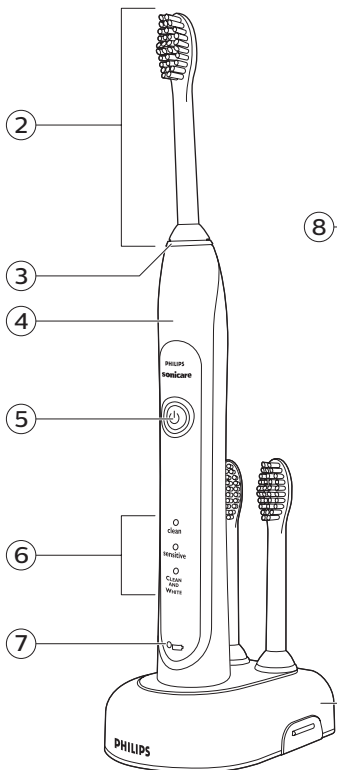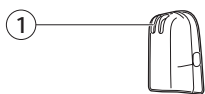
| | |
|-------------------------|-----------|
| English | 6 |
| Bahasa Indonesia | 18 |
| Bahasa Melayu | 30 |
| Français | 43 |
| Tiếng Việt | 57 |
| ภาษาไทย | 68 |
| فارسی | 89 |
| العربية | 100 |

Introduction

Congratulations on your purchase and welcome to Philips! To fully benefit from the support that Philips offers, register your product at www.philips.com/welcome.

IMPORTANT SAFEGUARDS

Warning: READ ALL INSTRUCTIONS BEFORE USE

DANGERS

To reduce the risk of electrocution:

- 1 Always place and store the charger at a safe distance from water.
- 2 Do not reach for a charger that has fallen into water. Unplug immediately.
- 3 Never use a charger with a damaged cord or plug.

WARNINGS

To reduce the risk of burns, electrocution, fire or physical injury:

- 1 Do not use attachments other than those recommended by the manufacturer.
- 2 Do not plug the charger into an outlet with a voltage other than specified on the bottom of the charger. This product is designed to operate within a range of 100 to 240 volts. Voltage converters do not warrant voltage compatibility.
- 3 Never force the plug into an outlet; if the plug does not easily fit into the outlet, discontinue use.
- 4 Discontinue use if the product appears damaged in any way (brush head, handle or charger).
- 5 The mains cord cannot be replaced. If the mains cord is damaged, discard the charger.
- 6 Keep the cord away from heated surfaces.
- 7 Do not use the charger outdoors.
- 8 This product contains no user-serviceable parts. Refer to 'Warranty and support' if the product no longer works properly or needs repair.

- 9 Use this product only for its intended use as described in this document or as recommended by your dental professional.
- 10 This product is not intended for use by persons (including children) with reduced physical, sensory, or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- 11 Children should be supervised to ensure that they do not play with the product.
- 12 Philips Sonicare is a personal care device and is not intended for use on multiple patients in a dental office or institution.

MEDICAL WARNINGS

- 1 Consult your dental professional before you use this product if you have had oral or gum surgery in the previous 2 months.
- 2 Contact your dental professional if excessive bleeding occurs after using this product or bleeding continues to occur after 1 week of use.
- 3 Discontinue use of this product and contact a physician/dentist if discomfort or pain is experienced.
- 4 This product complies with the safety standards for electromagnetic devices. If you have a pacemaker or other implanted device, contact your physician or the device manufacturer prior to use.
- 5 Consult your physician if you have other medical concerns.

Electromagnetic fields (EMF)

This Philips device complies with all applicable standards and regulations regarding exposure to electromagnetic fields.

SAVE THESE INSTRUCTIONS

Your Philips Sonicare (Fig. 1)

- 1 Hygienic travel cap
- 2 Brush head
- 3 Interchangeable colour code ring (for select models)

8 English

- 4 Handle with soft grip
- 5 Power on/off button
- 6 Brushing modes
- 7 Charge light indicator
 - Deluxe charge indicator (specific types only)
 - Standard charge indicator (specific types only)
- 8 Standard charger (specific types only)
 - a Charger cover with brush head holder
 - b Travel charger
- 9 Deluxe charger (specific types only)
 - a Charger cover with brush head holder
 - b Travel charger
 - c Charger base with cord wrap facility

Note: The content of the box may vary based on the model purchased.

Getting started

Changing the colour code ring (for select models)

Sonicare brush heads come with interchangeable colour code rings to identify your brush head.

To change the colour code ring:

- 1 Pull the colour code ring from the bottom of the brush head. (Fig. 2)
- 2 Install a new ring by slipping one edge of the new ring over the bottom of the brush head. (Fig. 3)

Attaching the brush head

- 1 Align the brush head so the bristles point in the same direction as the front of the handle. (Fig. 4)
- 2 Firmly press the brush head down onto the metal shaft until it stops.

Note: There is a small gap between the brush head and the handle.

Charging your Philips Sonicare

- 1 Plug the charger into an electrical outlet.
- 2 Place the handle on the charger. (Fig. 5)
 - The battery charge indicator next to the battery symbol flashes to indicate that the toothbrush is charging.

Battery charge indicator

The battery charge indicator shows the approximate remaining battery charge.

Deluxe charge indicator (specific types only):

- Solid green: 50-100%
- Flashing green: 10-49%
- Flashing yellow: less than 10%

Standard charge indicator (specific types only):

- Solid green: 50-100%
- Flashing green: less than 50%

Note: When the battery is low, you hear 3 beeps after the completion of the brushing cycle and the battery charge indicator flashes for 30 seconds.

Note: To keep the battery fully charged at all times, you may keep your Philips Sonicare on the charger when not in use. It takes at least 24 hours to fully charge the battery.

Using your Philips Sonicare

Brushing instructions

- 1 Wet the bristles and apply a small amount of toothpaste.
- 2 Place the toothbrush bristles against the teeth at a slight angle (45 degrees), pressing firmly to make the bristles reach the gumline or slightly beneath the gumline. (Fig. 6)

Note: Keep the center of the brush in contact with the teeth at all times.

- 3 Press the power on/off and mode button to turn on the Philips Sonicare.

10 English

- 4 Apply light pressure to maximize Philips Sonicare's effectiveness and let the Philips Sonicare toothbrush do the brushing for you. Do not scrub.
- 5 Gently move the brush head slowly across the teeth in a small back and forth motion so the longer bristles reach between your teeth. Continue this motion throughout the brushing cycle.
- 6 To clean the inside surfaces of the front teeth, tilt the brush handle semi-upright and make several vertical overlapping brushing strokes on each tooth.

Note: To make sure you brush evenly throughout the mouth, divide your mouth into 4 sections using the Quadpacer feature (see chapter 'Features').

- 7 Begin brushing in section 1 (outside upper teeth) and brush for 30 seconds before you move to section 2 (inside upper teeth). Continue brushing in section 3 (outside lower teeth) and brush for 30 seconds before you move to section 4 (inside lower teeth) (Fig. 7).

Note: When using the White mode, two additional 15 second intervals can be used for polishing your front teeth.

- 8 After you have completed the brushing cycle, you can spend additional time brushing the chewing surfaces of your teeth and areas where staining occurs. You can also brush your tongue, with the toothbrush turned on or off, as you prefer. (Fig. 8)

Your Philips Sonicare is safe to use on:

- Braces (brush heads wear out sooner when used on braces).
- Dental restorations (fillings, crowns, veneers).

Note: Ensure that braces or dental restorations are properly adhered and not compromised. If not, braces or dental restorations may be damaged when they come into contact with the brush head.

Note: When the Philips Sonicare is used in clinical studies, the default 2-minute Clean mode with high intensity must be selected. The handle should be fully charged and the Easy-start feature deactivated.

White mode brushing instructions

The Clean and White mode consists of 2 minutes of Clean mode to clean the whole mouth and an additional 30 seconds of White mode to focus on the visible front teeth.

- 1 Brush the first 2 minutes as explained in section 'Brushing instructions'.
- 2 After the 2 minutes of Clean mode, the White mode starts with a change in brushing sound and motion. This is your signal to start brushing the upper front teeth for 15 seconds.
- 3 At the next beep and pause, move to the bottom front teeth for the final 15 seconds of brushing (Fig. 9).

Brushing modes

Philips Sonicare automatically starts in the default Clean mode.

To select a different brushing mode:

- 1 Press the on/off button to toggle between modes.

The green LED indicates the selected mode.

Note: You can only switch between modes when the toothbrush is switched on.

- Clean mode: Standard mode for superior teeth cleaning.
- Sensitive mode (specific types only): Gentle yet thorough cleaning for sensitive gums and teeth.
- Clean and White mode: 2 minutes of Clean mode with an additional 30 seconds of White mode to focus on the visible front teeth.

Note: When the Sonicare is used in clinical studies, the default 2-minute Clean mode must be selected. However, in whitening/stain reduction studies, the Clean and White mode must be selected. The handle should be fully charged and the Easy-start feature deactivated.

Features

Easy-start

This Philips Sonicare comes with the Easy-start feature activated. The Easy-start feature gently increases the power over the first 14 brushings to help you get used to the brushing with the Philips Sonicare.

Note: Each of the first 14 brushings must be at least 1 minute in length to properly advance through the Easy-start ramp-up cycle.

Note: When the Clean and White mode is selected, the Easy-start feature is inactive.

Activating or deactivating the Easy-start feature

- 1 Attach the brush head to the handle.
- 2 Place the handle in the plugged-in charger.

To activate Easy-start:

Press and hold the power on/off and mode button for 2 seconds. You hear 2 beeps and the battery level indicator lights up green to indicate that the Easy-start feature has been activated.

To deactivate Easy-start:

Press and hold the power on/off and mode button for 2 seconds. You hear 1 beep and the battery level indicator lights up yellow to indicate that the Easy-start feature has been deactivated.

Note: Using the Easy-start feature beyond the initial ramp-up period is not recommended and reduces Philips Sonicare effectiveness in removing plaque.

Smartimer

The Smartimer indicates that your brushing cycle is complete by automatically turning off the toothbrush at the end of the brushing cycle. Dental professionals recommend brushing no less than 2 minutes twice a day.

If you want to pause or stop brushing during the 2-minute brushing cycle, press and hold the on/off button for 2 seconds. Press the on/off button again and the Smartimer picks up where you left off.

The Smartimer automatically resets to the beginning of a 2-minute cycle if:

- you stop brushing for 30 seconds or longer.
- you put the handle in the charger.

Quadpacer (select models)

- The Quadpacer is an interval timer that has a short beep and pause to remind you to brush the 4 sections of your mouth evenly and thoroughly. (Fig. 7)
- When you use the two-minute Clean or Sensitive mode, you hear a short beep and a pause in brushing action at 30, 60 and 90 seconds.
- When you use the Clean and White mode, you hear and feel a change in the brush head sound and motion after 2 minutes of Clean mode. Then move the brush head to the visible upper front teeth and brush them for 15 seconds. At the final beep and pause, move to the visible lower front teeth and brush them for 15 seconds. (Fig. 9)

Cord wrap facility (specific types only)

If your Sonicare type includes the deluxe charger cover and base, the travel charger is pre-installed inside the cover. If you want to shorten the mains cord, you can store excess cord in the cord wrap facility in the charger base.

- 1 To separate the charger cover from the charger base, press the two grey snap bars on the charger base and pull the white charger cover upwards. (Fig. 10)
- 2 Wrap excess mains cord around the grey charger base as shown in the picture. Make sure you wrap the cord on the inside of the two snap levers. (Fig. 11)
- 3 When you have stored the excess cord, guide the mains cord through the small groove in the rear of the grey charger base. (Fig. 12)
- 4 To reattach the charger cover, press it down over the charger base until it snaps into place.

Tip: For extra convenience during travel, you can remove the travel charger and use it without the charger cover and charger base.

Cleaning

Do not clean brush heads, the handle, the charger or the charger covers in the dishwasher.

Toothbrush handle

- 1 Remove the brush head and rinse the metal shaft area with warm water. Make sure you remove any residual toothpaste (Fig. 13).

Caution: Do not push on the rubber seal on the metal shaft with sharp objects, as this may cause damage.

- 2 Wipe the entire surface of the handle with a damp cloth.

Brush head

- 1 Rinse the brush head and bristles after each use (Fig. 14).
- 2 Remove the brush head from the handle and rinse the brush head connection with warm water at least once a week.

Charger

- 1 Unplug the charger before you clean it.
- 2 Wipe the surface of the charger with a damp cloth. Make sure you remove any toothpaste or other residue from the charger surfaces.

Storage

- If you are not going to use the Sonicare for an extended period of time, unplug the charger from the electrical outlet, clean it and store it in a cool and dry place away from direct sunlight. You can use the cord wrap facility (specific types only) to store the mains cord in an organised way.
- You can store brush heads on the pegs at the back of the charger.

Replacement

Brush head

Replace Philips Sonicare brush heads every 3 months to achieve optimal results. Use only Philips Sonicare replacement brush heads.

Recycling

- Do not throw away the product with the normal household waste at the end of its life, but hand it in at an official collection point for recycling. By doing this, you help to preserve the environment.

Removing the rechargeable battery

Warning: Only remove the rechargeable battery when you discard the appliance. Make sure the battery is completely empty when you remove it.

To remove the rechargeable battery, you need a flat-head (standard) screwdriver. Observe basic safety precautions when you follow the procedure outlined below. Be sure to protect your eyes, hands, fingers, and the surface on which you work.

- 1 To deplete the rechargeable battery of any charge, remove the handle from the charger, turn on the Philips Sonicare and let it run until it stops. Repeat this step until you can no longer turn on the Philips Sonicare.
- 2 Insert a screwdriver into the slot located in the bottom of the handle and turn counterclockwise until you see a gap between the bottom cap and the handle. (Fig. 16)
- 3 Insert the screwdriver into the gap and pry the bottom cap out of the handle.
- 4 Hold the handle upside down and push down on the shaft to release the internal components of the handle. (Fig. 17)
- 5 Insert the screwdriver under the circuit board, next to the battery connections, and twist to break the connections. Remove the circuit board and pry the battery from the plastic carrier. (Fig. 18)

The rechargeable battery can now be recycled or disposed of and the rest of the product discarded appropriately.

Guarantee and support

If you need information or support, please visit www.philips.com/support or read the separate worldwide guarantee leaflet.

Guarantee restrictions

The terms of the international guarantee do not cover the following:

- Brush heads.
- Damage caused by use of unauthorised replacement parts.
- Damage caused by misuse, abuse, neglect, alterations or unauthorised repair.
- Normal wear and tear, including chips, scratches, abrasions, discolouration or fading.

Frequently asked questions

This chapter lists the questions most frequently asked about the appliance. If you cannot find the answer to your question, contact the Consumer Care Center in your country.

| Question | Answer |
|---|--|
| Why doesn't the Sonicare work? | You may need to recharge the Sonicare. Charge the Sonicare for at least 24 hours. |
| | The socket may not be live. The electricity supply to sockets in mirror cabinets over washbasins or shaver sockets in bathrooms may be cut off when the light is switched off. |
| | You may need to reattach or replace the brush head. |
| Why do I feel a tickling sensation when I use the Sonicare? | You may experience a slight tickling/tingling sensation when you use the Sonicare for the first time. As you become accustomed to brushing with the Sonicare, this sensation diminishes. |
| The brush head chatters against my teeth. What can I do about this? | To avoid this, keep the brush head properly oriented, with the bristles placed on the gumline at a slight angle. |
| Why does the Sonicare seem less powerful? | You may need to charge the Sonicare. |

| Question | Answer |
|---|---|
| | The Easy-start feature may be on. Deactivate the Easy-start feature (See chapter 'Features'). |
| | You may need to replace the brush head. |
| Why does the Sonicare seem to change power during brushing? | You may be in the Clean and White brushing mode. If this is the case, the toothbrush starts in the Clean mode for 2 minutes and then changes to the White mode for 30 seconds. The Clean and White mode has a different sound and brush movement. |
| Do I need a special charger if I travel with the Sonicare? | No, you do not need a special charger. The Sonicare charger can be used with common worldwide voltages between 100–240V AC, at either 50 or 60Hz. A simple plug adapter is all that is required. |

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©2016 Koninklijke Philips N.V. (KPNV). All rights reserved.
Philips and the Philips shield are trademarks of KPNV.
Sonicare and the Sonicare logo are trademarks
of Philips Oral Healthcare, LLC and/or KPNV.
Healthcare, LLC et (ou) KPNV.