
**sonicare**<sup>®</sup>  
the **sonic** toothbrush

elite  
e5000  
series

**PHILIPS**
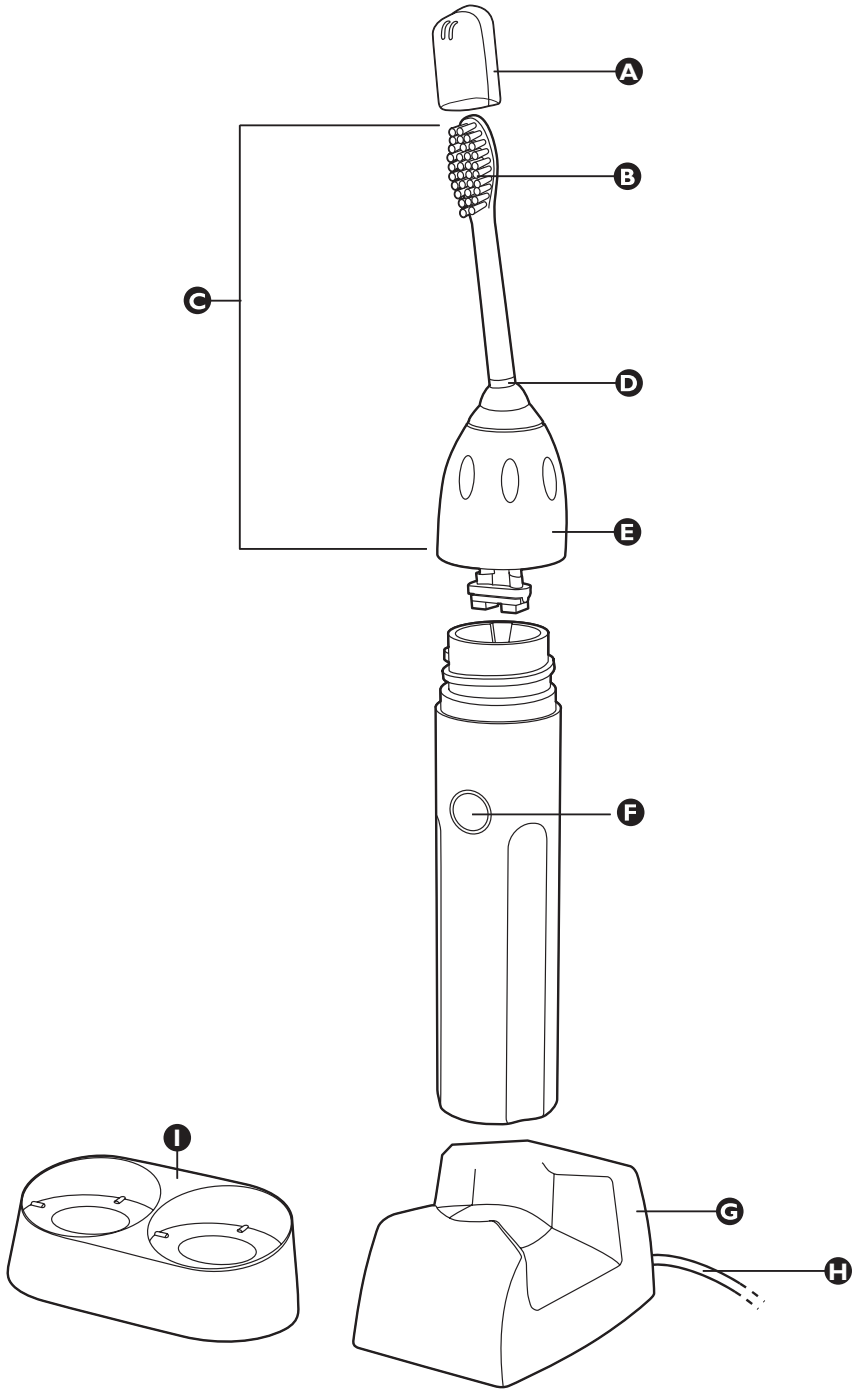
<b>ENGLISH</b>	<b>4</b>
<b>DANSK</b>	<b>14</b>
<b>ESPAÑOL</b>	<b>24</b>
<b>SUOMI</b>	<b>34</b>
<b>NORSK</b>	<b>44</b>
<b>PORTUGUÊS</b>	<b>54</b>
<b>РУССКИЙ</b>	<b>64</b>
<b>SVENSKA</b>	<b>74</b>

**General description (Fig. 1)**

- A** Travel cap
- B** Extra-soft bristles
- C** Brush head
- D** Colour code ring
- E** Removable nut
- F** On/off button & charging light
- G** Charger with cord storage
- H** Mains cord
- I** Brush head holder


**Important**

Read these instructions for use carefully before you use the appliance and save them for future reference.

- ⓘ Check if the voltage that is indicated on the bottom or back of the charger corresponds to the local mains voltage before you connect the charger.
- ⓘ Voltage converters and plug adapters **DO NOT** guarantee voltage compatibility. Always check for voltage compatibility.
- ⓘ Keep the charger away from water. Do not place or store it over or near water that is contained in a bathtub, washbasin, sink etc. Do not immerse the charger in water or any other liquid.
- ⓘ Never use the charger if it is damaged or if the cord or plug is damaged. If the charger, the cord or the plug is damaged, always have the charger replaced by one of the original type in order to avoid a hazard.
- ⓘ Unplug the charger before you clean it and make sure that it is dry before you put the plug in the wall socket.
- ⓘ Only use the charger supplied.
- ⓘ Do not use the charger outdoors.
- ⓘ Do not use the toothbrush in the bath or shower.
- ⓘ Close supervision is necessary when this appliance is used by or on children under the age of 10 or the disabled.
- ⓘ Stop using the appliance if you are unable to handle it properly.
- ⓘ If your toothpaste includes peroxide, baking soda or other bicarbonate (common in whitening toothpaste), thoroughly clean the brush head with soap and water after each use.
- ⓘ If you have had oral or gum surgery in the previous 2 months, consult your dentist before you use the toothbrush.
- ⓘ Consult your dentist if excessive bleeding occurs after using this toothbrush or if bleeding continues to occur after 1 week of use.
- ⓘ The Sonicare toothbrush complies with the safety standards for electromagnetic devices. If you have a pacemaker or other implanted device, consult your doctor or the manufacturer of the implanted device if you have any questions about using the Sonicare.
- ⓘ Do not use other brush heads than the ones the manufacturer recommends.
- ⓘ Replace the brush head every 6 months or sooner if signs of wear appear, such as crushed or bent bristles.
- ⓘ This appliance has only been designed for cleaning teeth, gums and tongue. Do not use it for any other purpose.


- ▮ Supervise children when they are brushing until they have established good brushing habits.
- ▮ The Sonicare toothbrush is a personal care device and it is not intended for use on multiple patients in a dental practice or institution.
- ▮ The mains cord of the charger cannot be replaced. If the mains cord is damaged, you must discard the charger.

### Preparing for use

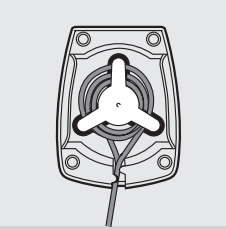
- 1 Store excess cord by winding it round the cord storage hooks in the bottom of the charger.

### Attaching the brush head

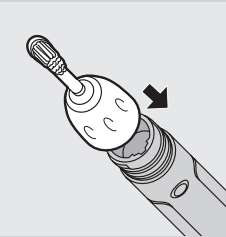
- 1 Place the brush head onto the handle by sliding the ribs on the inside of the nut into the grooves in the top of the handle.

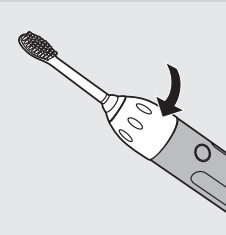- 2 Screw the brush head nut clockwise until it is firmly tightened. If the brush head rattles during use, tighten the nut further.
- 3 Remove the travel cap from the brush head before you start to brush.


### Charging

Charge the Sonicare for at least 24 hours during the first few days of use. We advise you to keep the Sonicare in the plugged-in charger when not in use to maintain a full battery charge. When the battery is empty, it takes at least 24 hours to fully charge the Sonicare.

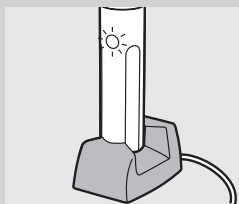- 1** To charge the Sonicare, place the handle properly in the charger with the on/off button facing towards the front of the charger.

Make sure the bottom of the handle sits back securely in the charger so that it makes complete contact.

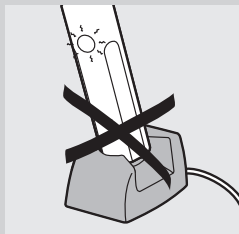- When the handle is placed incorrectly in the charger, the charging light may still go on, but charging does not take place.

The handles of other Sonicare models do not fit in this charger and therefore will not charge effectively. Do not attempt to force a handle into this charger, as this will cause damage.

## Using the appliance


### Useful tips

- The brush head comes with a convenient travel cap. Remove the travel cap from the brush head before you switch on the Sonicare.
- Place the brush head in your mouth before you switch on the appliance and keep your mouth closed during brushing to prevent splattering of toothpaste.
- A fully charged Sonicare allows you to brush your teeth for 2 minutes twice a day for approximately 2 weeks without charging. This makes it easy to travel without the charger.
- Voltage converters and plug adapters **DO NOT** guarantee voltage compatibility. Always check for voltage compatibility.

### Brushing technique

Follow these simple steps to maximise your brushing experience.

- 1** Before you switch on the Sonicare, place the bristles of the brush head on your gumline at a 45-degree angle.
- 2** Hold the handle with a light grip and switch on the Sonicare. Gently move the bristles in a small circular motion so that the longer bristles reach between your teeth. After a few seconds, pass the bristles to the next section. Continue this routine throughout the 2-minute brushing cycle.
- 3** To maximise the effectiveness of the Sonicare toothbrush, use light pressure and let the Sonicare do the brushing for you. Do not scrub.

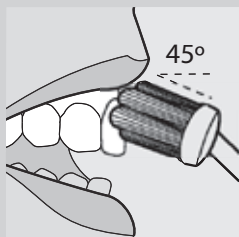- 4** To make sure you brush evenly throughout the mouth, divide your mouth into 4 sections: outside upper teeth, inside upper teeth, outside lower teeth and inside lower teeth. Begin brushing the first section (outside upper teeth) and brush for 30 seconds before moving to the next section (inside upper teeth). Continue brushing the other two sections for 30 seconds each until you have completed the 2-minute brushing cycle.

- 5** After the 2-minute brushing cycle, you can briefly switch on the Sonicare again to brush the chewing surfaces of your teeth.

Press the on/off button to switch the Sonicare on or off as needed.

To clean your mouth more completely, you can also brush your tongue with the toothbrush switched on or off, as you prefer.

### Tips for specific oral care needs

#### ► Stain removal

For even better stain removal, dedicate a few seconds of extra brushing time to areas where excess staining occurs.

#### ► Braces

Place the brush head at a 45-degree angle between the brackets and the gumline. Using light pressure, gently move the bristles in a small circular motion so the longer bristles reach between your teeth for 2 seconds.

Then gently roll the brush handle so the bristles clean the bracket and the rest of the tooth. Subsequently pass the bristles to the next tooth.


Continue this procedure until you have brushed all your teeth. Please note that brush heads wear out sooner when used on braces.

#### ► Dental restorations

Use the Sonicare on dental restorations (such as implants, caps, crowns, bridges and fillings) as you do on your natural tooth surfaces. Work the bristles round all the areas of the restoration for thorough cleaning.

#### ► Periodontal pockets

Your dentist or dental hygienist may have indicated that you have periodontal pockets, deep pockets where plaque bacteria collect. If this is the case, spend extra time brushing these areas after you have completed the 2-minute brushing cycle to improve your gum health.
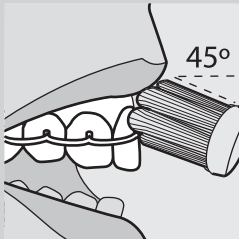
## Features

### Easy-start

- All Sonicare models come with the Easy-start feature activated, which helps you get used to brushing with the Sonicare.
- The brushing power gently increases to full power over the first 12 brushings.
- The Easy-start feature is active if you hear 2 beeps immediately after you switch on the Sonicare. Once the Sonicare has reached optimal brushing power (after the 12th brushing), you no longer hear the 2 beeps at the beginning of each brushing cycle.

*Note: Each of the first 12 brushings must be at least 1 minute long to move properly through the Easy-start ramp-up cycle.*

To deactivate or activate the Easy-start feature:

- 1** Attach the brush head to the handle.

- 2 Place the handle in the plugged-in charger.
- 3 To deactivate the Easy-start feature, press and hold the on/off button for 5 seconds. You hear 1 beep to indicate that the Easy-start feature has been deactivated.
- 4 To activate the Easy-start feature, press and hold the on/off button for 5 seconds. You hear 2 beeps to indicate that the Easy-start feature has been activated.


Use of the Easy-start feature beyond the initial ramp-up period is not recommended and reduces the effectiveness of the Sonicare toothbrush in removing plaque.

Make sure the Easy-start feature is deactivated and the handle fully charged when the Sonicare is used in clinical trials.

### Smartimer

- ▶ The Smartimer automatically sets the brushing time to the 2 minutes that dentists and dental hygienists recommend.
- ▶ After 2 minutes, the Smartimer automatically switches off the toothbrush.
- ▶ If you want to stop brushing briefly during a 2-minute brushing cycle, press the on/off button. When you press the on/off button again, the Smartimer resumes the brushing cycle at the point where you stopped.
- ▶ The Smartimer automatically resets to the beginning of a 2-minute cycle if:
  - you stop brushing for 30 seconds or longer;
  - you put the handle in the charger.

### Quadpacer® (specific models only)

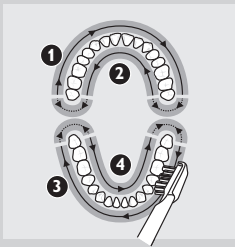
- ▶ The Quadpacer 30-second interval timer ensures that you brush all sections of your mouth evenly and thoroughly. At 30, 60 and 90 seconds, you hear a short beep and notice a pause in brushing action. These signals tell you to continue with the next section of your mouth.

### Battery charge indication (specific models only)

- ▶ When the battery is running low, you hear 5 beeps after the 2-minute brushing cycle and the charging light blinks rapidly. When you hear the low charge indication, there are 3 or fewer 2-minute brushings left.

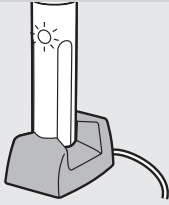- 1 Place the handle in the charger (with the on/off button facing towards the front of the charger) to recharge the appliance.
- ▶ The charging light blinks slowly when the handle is in the charger to indicate that charging is in progress.
- ▶ When the Sonicare is fully charged, the charging light burns continuously to indicate that the Sonicare is ready for use.

Keep the Sonicare in the charger when you are not using it to maintain a full battery charge.

## Cleaning

Regular cleaning guarantees the best performance and a long life for the appliance.

- 1 Rinse the brush head and bristles after each brushing.

- 2 Clean the brush head at least once a week, including under the removable nut.


To remove the nut:

- ▶ Gently press on the rubber seal with your thumbs.

Do not use sharp objects to press on the seal, as this may cause damage.

- ▶ Pull the nut up over the top of the brush head.

To reattach the nut, slip the nut over the brush head shaft and pull it down until it snaps home.

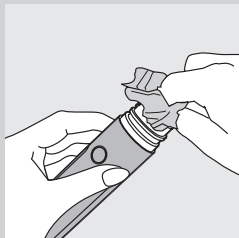**3** Clean and dry the handle weekly.

**4** Clean the charger periodically.

Unplug the charger before you clean it. After you have cleaned the charger, dry it completely before you put the plug back into the wall socket.

We advise you to use a mild cleaning agent, water and a moist cloth to clean the Sonicare.

### Storage

If you do not intend to use the appliance for a long time, unplug it, clean it and store it in a cool and dry place away from direct sunlight. You can store the cord by winding it round the cord storage hooks in the bottom of the charger.

### Replacement

#### Brush head

- ▶ Replace the brush head every 6 months to achieve optimal results.
- ▶ Only use Sonicare replacement brush heads that are suitable for this Sonicare model.

### Environment


- ▶ Do not throw away the appliance with the normal household waste at the end of its life, but hand it in at an official collection point for recycling. By doing this, you help to preserve the environment.
- ▶ The built-in rechargeable batteries contain substances that may pollute the environment. Always remove the batteries before you discard and hand in the appliance at an official collection point. Dispose of the batteries at an official collection point for batteries. If you have trouble removing the batteries, you can also take the appliance to a Philips service centre. The staff of this centre will remove the batteries for you and will dispose of them in an environmentally safe way.

### Removing the batteries


This process is not reversible.

To remove the batteries yourself, you need a flat-head (standard) screwdriver:


- 1** Remove the toothbrush from the charger, switch it on and let it run until it stops. Continue to switch on the Sonicare until the battery is completely empty.

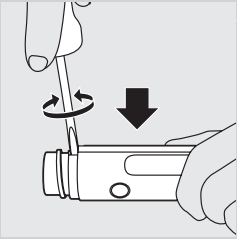
- 2** Insert a screwdriver into one of the slots that are located near the top part of the handle. Push the screwdriver blade into the handle as far as it goes and twist the blade. The housing begins to split.

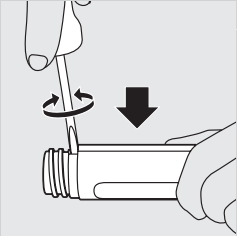
- 3** Remove the screwdriver and apply the above instruction to the other slot, which is located on the opposite side of the handle. Twist the blade to split the rest of the housing.

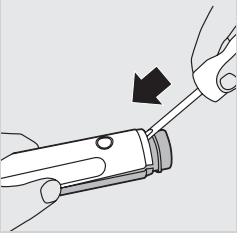- 4** Now insert the screwdriver into the middle slot, which is located on the top front part of the handle, as shown. Push the blade in as far as it goes.

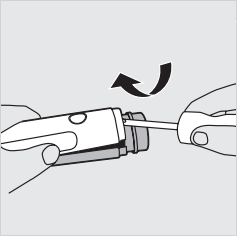
- 5** Now use the screwdriver as a lever to pry the 2 handle sections apart.

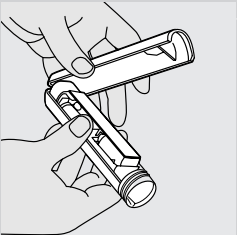
- 6** Pull the 2 handle sections completely apart.


- 7** Insert the screwdriver underneath the electronic circuit board in multiple places.

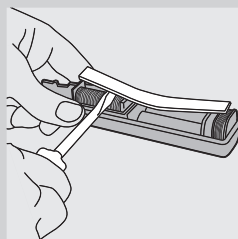- 8** Twist and pry the electronic circuit board loose until it is completely separated from the handle.

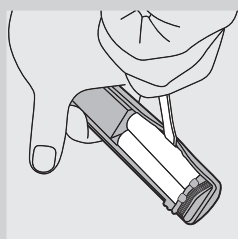- 9** Pry the batteries loose by inserting the screwdriver on the sides of the two batteries in the locations that are shown in the figure.

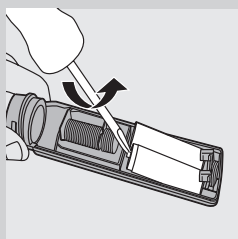- 10** Twist the screwdriver to break the glue bond round the batteries. The 2 batteries are joined together with glue and will be lifted out together.

### Guarantee & service

If you need information or if you have a problem, please visit the Philips website at [www.philips.com](http://www.philips.com) or contact the Philips Customer Care Centre in your country (you find its phone number in the worldwide guarantee leaflet). If there is no Customer Care Centre in your country, turn to your local Philips dealer or contact the Service Department of Philips Domestic Appliances and Personal Care BV.

### Guarantee restrictions

The terms of the international guarantee do not cover the following:

- Brush heads
- Damage caused by misuse, abuse, neglect or alterations

- Normal wear and tear, including chips, scratches, abrasions, discolouration or fading

**Troubleshooting**

This chapter summarises the most common problems you could encounter with the Sonicare. If you are unable to solve the problem, contact the Philips Customer Care Centre in your country (see chapter 'Guarantee & service') or visit [www.philips.com](http://www.philips.com).

Problem	Solution
The Sonicare toothbrush does not work.	You may need to recharge the Sonicare. Charge the Sonicare for at least 24 hours.
	Make sure the socket is live when you recharge the appliance. The electricity supply to sockets in mirror cabinets over washbasins or shaver sockets in bathrooms may be cut off when the light is switched off.
The motor runs but the brush head does not operate.	You may need to reattach or replace the brush head.
I feel a tickling sensation when I use the Sonicare.	You may experience a slight tickling/tingling sensation when you use the Sonicare for the first time. As you become accustomed to brushing with the Sonicare, this sensation diminishes.
The brush head rattles when I use the Sonicare.	The Sonicare should produce a gentle hum during operation. Tighten the brush head nut to eliminate excessive rattling. If the rattling does not stop, you may need to replace the brush head.
	Make sure you remove the travel cap before you switch on the appliance.
The brush head chatters against my teeth.	You can avoid this by keeping the brush head properly oriented, with the brush head placed on the gumline at an angle of 45 degrees.
Brushing seems less powerful.	You may need to recharge the Sonicare.
	The Easy-start feature may be active and needs to be deactivated (see chapter 'Features').
	You may need to clean the brush head.
	You may need to replace the brush head.
The Sonicare gives off an unpleasant smell.	You may need to clean the appliance. Clean the handle, brush head and charger thoroughly (see chapter 'Cleaning').
	You may need to replace the brush head.

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© 2006 Philips Oral Healthcare, Inc. All rights reserved. PHILIPS and the Philips shield are registered trademarks of Koninklijke Philips Electronics N.V. Quadpacer, Sonicare and the Sonicare logo are trademarks of Philips Oral Healthcare, Inc. and/or KPENV.
Philips Oral Healthcare, Inc.
U.S.A.
Printed in U.S.A.