

# PHILIPS

IPL - Hair removal device

Lumea Prestige

Prevent hair regrowth

With SmartSkin sensor

3 intelligent attachments

Corded and cordless use

BG9041/83



## Prevent hair regrowth\*

with Intense Pulsed Light (IPL) technology

Philips Lumea uses Intense Pulse Light (IPL) technology for safe, long-term hair free. This effective professional salon technology prevents hair regrowth at the root, suitable for home use.

### Continuous smooth skin with Lumea IPL technology

- IPL technology to prevent hair regrowth at the root
- Suitable hair and skin types
- Two phases treatment: initial phase and touch-up phase
- Two flash modes: slide and stamp

### Maximum freedom, minimum hassle

- 5 adjustable light intensity settings
- Corded and cordless use

### Uniquely designed for specific body areas

- Uniquely curved attachments adapts to each body area
- Body attachment with curved-in design for large body areas
- Underarm attachment with curved-out design
- Sensitive area attachment with special designed extra filter

# Highlights

## Long lasting hair removal



Philips Lumea uses an innovative light based technology, called IPL, to break the cycle of hair regrowth. Gentle pulses of light are applied at the root, putting the follicle into a resting phase. As a consequence, the amount of hair your body grows gradually decreases. Developed closely with skin experts and dermatologists, Philips has adapted this technology for safe and effective use in the comfort of your home.

## SmartSkin sensor



Philips Lumea Prestige has 5 energy settings and makes it easy for you to choose one thanks to the SmartSkin sensor.


## Suitable hair and skin types



Philips Lumea works effectively on a wide variety of hair and skin types. Works on (naturally) dark blonde, brown and black hairs and on skin tones ranging from very white (I) to dark brown (V). As with other light based


treatments, Philips Lumea is NOT meant to treat red, light-blonde or white/grey hair. This product is also NOT suitable for very dark skin (skin type VI).

## Treatment schedule

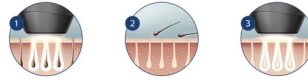

For best results, we recommend an 8-week schedule, followed by touch ups. 1. Initial phase: apply Philips Lumea every 2 weeks for the initial 4 – 5 treatments; 2. Touch-up phase: afterwards, repeat the treatment every 4 – 8 weeks.

## Application body areas

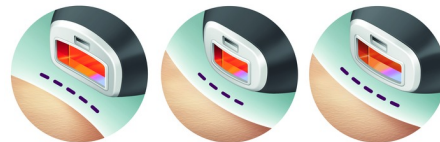

Philips Lumea has been developed with dermatologists for a safe, effective and gentle application on arms, underarms, shoulders, stomach, pubic area (excluding the genital area and scrotum), sensitive areas, legs, chest and back.

## Body attachment

For body area — 4.1cm<sup>2</sup> window size. The body attachment is curved in for optimal skin contact when treating large areas like legs.

## Underarm attachment




For underarm area — 3cm<sup>2</sup> window size. The specially curved attachment is designed to treat hard-to-reach underarm hair.

## Sensitive area attachment




For sensitive area — 3cm<sup>2</sup> window size. The curved attachment has an extra transparent filter and specially designed for sensitive area.

## Corded and cordless use

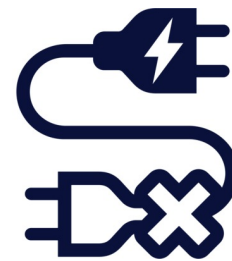

This Philips Lumea is the only IPL device that can be used with or without cord for complete convenience. Use cordless for full freedom of movement and plug in when needed for treatment with no interruptions.

## Slide and flash mode

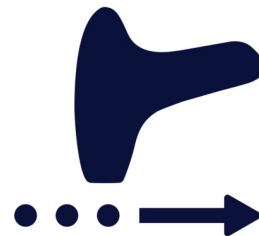

The slide and flash modes enables fast treatment on different body areas.

# Specifications

## Technical specs. attachments

**Body treatment attachment:** Window size: 4.1 cm<sup>2</sup>, Shape: Convex curved, Tailored treatment for body

**Underarm treatment attachment:** Window size: 3cm<sup>2</sup>, Shape: Concave curved, Tailored for underarm area

**Sensitive treatment attachment:** Window size: 3cm<sup>2</sup>, Shape: Concave curved, extra filter, Tailored for sensitive area

## Technical specifications

**High performance lamp:** Delivers > 250,000 flashes

## Application areas

Body areas

Sensitive areas

Underarm

## Items included

**Instructions for use:** User manual

**Storage:** Luxurious pouch

**Adapter:** 19.5V / 4000mA

## Power

**Voltage:** 100-240 V

## Safety and adjustable settings

**Integrated UV filter:** Protects skin from UV light

**5 light energy settings:** Adjustable to your skin type

**Integrated safety system:** Prevent unintentional flashing

**Skin Tone sensor:** Detects your skin tone

**SmartSkin sensor:** The right setting on demand

## Service

**Warranty:** 2 years global warranty

## Application mode

**Slide and Flash:** For quick application

**Corded / cordless use:** Both cordless and corded use

**Stamp and Flash:** For treatment on small areas

