



Pediatric compressor nebulizer system

- Home use table-top device
- Treatment in 6-8 minutes*
- *Using a 2.5ml Salbutamol dose
- Designed for children

HH1304/00

Making therapy fun

Sami the Seal pediatric nebulizer compressor system is shaped like a friendly animal character to help encourage your child to take their nebulizer therapy.

Benefits

Specially designed for children

- Efficient, reliable and friendly
- Cute and effective

Short treatment times

- Treatment in 6-8 minutes*

Proven SideStream technology

- Consistently produces high-quality aerosol
- Choose which nebulizer is right for you and your child

For commonly prescribed medications

- Compatible with commonly prescribed medications

Features

Making therapy fun




The Sami the Seal compressor nebulizer system combines power and reliability with a child-friendly design and cute character to provide you and

your child with a more user- friendly treatment experience. The compressor is bottom-weighted to provide extra stability.

PHILIPS

Sami the Seal

Tucker pediatric mask

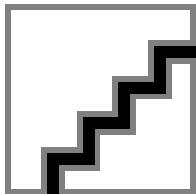

Our cute Tucker the Turtle character mask is made of a soft, flexible material that contours to your child's face for comfort and fit. It also incorporates a vent design that helps to direct medication away from your child's eyes.

Fast, effective treatment



Our SideStream nebulizer is specially designed to keep your child's treatment times as short as possible. As Sami the Seal pushes air through the tube, additional air is drawn in through the top of the nebulizer to create a fast flow of aerosol medication for your child to inhale. Treatments can take as little as 6 minutes*

Proven nebulizer technology




Consistent delivery of high quality aerosol therapy. SideStream is designed to withstand frequent, repeated use and is able to efficiently deliver your child's nebulizer medication until it is time for replacement.

Multiple configurations



The system is available in two configurations – one for use with the SideStream disposable nebulizer only, and one with the SideStream disposable and SideStream reusable nebulizers.

Versatility

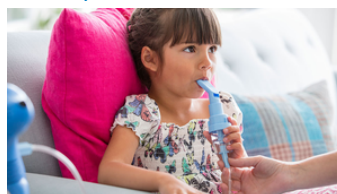

Sami the Seal is designed to use with most aerosol medications prescribed for common respiratory diseases.

Specifications

Please note that this is a pre-sales leaflet. The contents of this leaflet reflect the best of our knowledge per date and country mentioned above. The contents of this leaflet are subject to change without notice. Philips does not accept any liability as to the contents of this leaflet.

Design

Reliable and appealing Child-friendly

Operating parameters

Maximum pressure	29-43.5 psi
Operating mode	Continuous
Performance average flow rate	6 LPM @ 10 psig
Thermal fuse	Functioning temp 130 degrees C

Product details

Weight	3.5 lbs
Size	226mm x 155mm x 267mm
Warranty	5 years
Treatment time (6-8 minutes*)	* 2.5ml Salbutamol dose
Electrical requirements	120VAC-60 Hz
Power	120 VA

** Using a 2.5ml Salbutamol dose



data subject to change
2018, January 20

Version: 11.1.1
UPC: 383730932682

© 2018 Koninklijke Philips N.V.
All rights reserved.

Specifications are subject to change without notice. Trademarks are the property of Koninklijke Philips N.V. or their respective owners.

www.philips.com