

# PHILIPS

Minimal contact  
nasal mask

DreamWear

Pack

and 4 cushion sizes

3 frame sizes

HH1121



## Like wearing nothing at all\*

With the freedom and comfort of DreamWear, you can sleep how you want \*2015 Philips User Preference Questionnaire

### Open field of vision

- Open design offers an extremely wide field of vision

### Innovative design

- Air tube connection at the top of the head
- Under-the-nose cushion
- Supports your sleep therapy
- Four-part modular design

### Comfortable design

- Slim-line headgear
- Soft fabric wraps
- Soft, flexible frame

# Highlights

## Freedom of movement



By moving the tube location to the top of the head, DreamWear offers more freedom of movement during the night. Users stated they didn't need to choose a sleep position when wearing DreamWear.\*

## Four-part modular design



Modular design for quick assembly and cleaning, with interchangeable cushions available in four sizes.

## Slim-line headgear




DreamWear combines comfort with aesthetics: Not only will you barely feel the mask, it is also subtle and unobtrusive.

## Soft fabric wraps



Provides a soft interface to the cheeks

## Soft, flexible frame

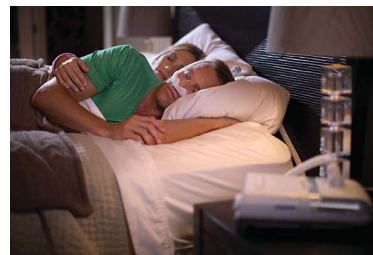

The flexible material through which the air flows is applied gently to the skin and is pleasant to wear.

## Supports your sleep therapy




DreamWear combines the advantages of minimal contact masks and nasal masks. It is comfortable and can help a successful start to sleep therapy.

## Under-the-nose cushion



Prevents red marks, discomfort or irritation in the nostrils or on the nose bridge. Users rated DreamWear more comfortable than their prescribed mask.

## Unobstructed view

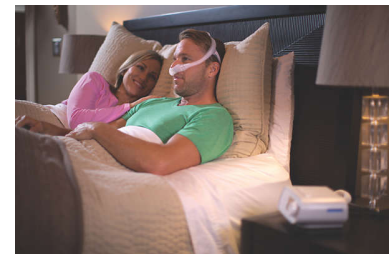

It's easy to wear glasses, read, watch TV, and use a computer or tablet before falling asleep.

# Specifications

## Innovative design

**Open field of vision:** Wear glasses, read or watch TV

## Specifications

**Operating pressure:** 4 -20 cmH<sub>2</sub>O

**Sound pressure level:** 18 dBA

**N° of parts:** Five

**Cushion material:** Silicone

**Headgear material:** Urethane foam/nylon/Spandex

**Frame material:** Polycarbonate

**Cleaning instructions: mask:** Hand wash with mild detergent, Rinse thoroughly, air dry

**Cleaning: headgear, tubing:** Hand wash with mild detergent, Rinse thoroughly, hang to dry

**Frame/Elbow/Swivel material:** Polycarbonate

**Fabric sleeves material:** Polyester/Spandex

## Comfort

**Under-the-nose cushion:** Prevents red marks in nostrils

## Intent of use

**To be used on patients <gt/>30kg:** with CPAP prescription

## Service

**Warranty:** 90-days

