



Philips DreamStation Go Travel CPAP machine

Same PAP performance on the go

Easy to use

Weighing less than 2 pounds
12mm micro-flexible tubing

HH1511



Reimagining the travel CPAP machine

Our sleek design has 1/3 fewer components to pack than other travel PAPs,* and DreamStation Go offers the same therapy and comfort features enjoyed by over 5 million PAP users.

Philips smallest PAP machine

- Built-in USB charging port
- Compact footprint
- Quick, simple setup and easy navigation

Philips lightest PAP at less than two pounds

- Overnight battery (Sold separately)

Philips smallest and lightest tubing**

- Micro-flexible tubing

Complies to FAA requirements*

- For airplane travel*

Store data for your physicians and clinicians

- Bluetooth is standard with every device
- Compatibility with Philips DreamMapper
- MicroSD card sold separately.

PHILIPS

Travel CPAP machine

Same PAP performance on the go Easy to use, Weighing less than 2 pounds, 12mm micro-flexible tubing

HH1511/00


Highlights

Built-in power supply



Built in power supply reduces the number of parts to carry. With these dimensions of 15 cm L X 15 cm W x 5.8 cm H, it fits easily even into a small bag, light and discreet. It can be carried anywhere. A pouch for mask and tubing is included in the kit.

Built-in USB charging port

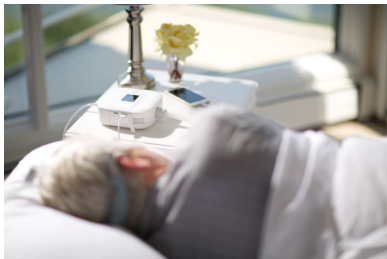

Allows for easy charging of cell phones or other devices, leaving outlets available for lamps, alarm clocks, etc. Convenient for travel.

Complies to FAA requirements



For your convenience at security stations in airports, there is a note on the bottom of the device stating that it is medical equipment and is suitable for airline use. The device is suitable

for use on airlines when it is operating from an AC power source or battery pack. Contact the airline prior to departure for information on power source available in flight. This device automatically compensates for altitude up to 7,500 feet. No manual adjustment is necessary.

Intuitive color touchscreen



Adjust device settings and view therapy information with DreamStation Go's colourful and intuitive touch screen. Take advantage of our MyStart feature that allows you to set the starting pressure for Auto CPAP therapy that is most comfortable for you. If the default air pressure feels too high or too low while trying to fall asleep you can simply adjust the MyStart pressure setting up or down, allowing you to fall asleep more comfortably.

Micro-flexible tubing



Included with your DreamStation Go purchase, the new 12mm micro-flexible tubing is Philips smallest and lightest we've ever offered, providing enhanced packability, flexibility and freedom of movement while you sleep. Fast and easy click connections on each

end of the tubing. Features a slim connection to our DreamWear mask, and includes a 22mm ISO adapter for other masks.

Overnight battery



If you plan on traveling off the power grid, DreamStation Go's overnight battery can help you stay committed to your therapy. DreamStation Go set at 10cm of pressure displayed an average run time of 13 hours*** when running off the Dreamstation Go overnight battery. Battery run times may vary due to higher CPAP pressure settings, excessive mask leakage, altitude, temperature, charging devices with USB port and others.

Bluetooth



Allows the ability to track and monitor your combined progress across both home and portable PAPs.

Specifications

Service

- Warranty: 2 years

Technical information

- Pressure range: 4 to 20 cm H₂O
- Pressure mode(s): CPAP, APAP
- Flex pressure relief: 1 to 3
- Flex modes: C-Flex, C-Flex+, A-Flex
- Ramp time: 0 to 45 min (5-min increments)
- Starting ramp pressure: Patient adjustable, Smart Ramp
- Dimensions: 15 cm L x 15 cm W x 5.8 cm H
- Weight: 844 g (1.86 lb)
- Filters: Reusable pollen
- Device setup: Touchscreen color display
- Data storage capacity: MicroSD card (optional), On-board: > 1 year
- Compliance meter: Breathing detection
- Altitude compensation: Automatic
- Electrical requirements: 100-240 VAC, 50/60 Hz



Issue date 2018-01-10

Version: 1.1.1

12 NC: 000B LG5 02S15
UPC: 6 06959 04846 1

© 2018 Koninklijke Philips N.V.
All Rights reserved.

Specifications are subject to change without notice.
Trademarks are the property of Koninklijke Philips N.V.
or their respective owners.

www.philips.com

* When compared against HDM Z1 and Somnetics Transcend devices (2016).

** When compared to Philips 22mm Performance tubing.

***PAP run-time test conditions: CPAP mode, pressure 10 cm H₂O, 12mm tubing, 37LPM of leak, 73.4 F room at an elevation of 1650 ft.