




Philips DreamWear
Nasal mask

No red marks on bridge of nose

Under-the-nose nasal cushion
Therapy device sold separately

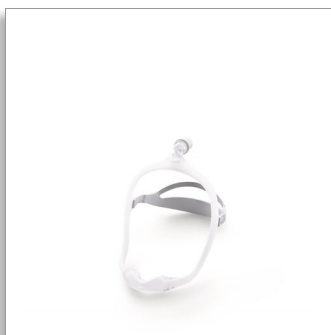

HH1135/00

Like wearing nothing at all*

The DreamWear All in One Starter is designed with freedom and comfort in mind. The pack comes with 3 frame and 4 cushion sizes. If you are not sure of your size, the DreamWear Nasal All in One Starter Pack is right for you.

Magnetic clips

- The headgear clips and mask cushion contain magnets

Open field of vision

- Open design offers an extremely wide field of vision

Innovative under-the-nose design

- Comfort and performance

Innovative design

- Air tube connection at the top of the head
- Four-part modular design

Comfortable design

- Soft fabric wraps
- Soft, flexible frame

PHILIPS

Nasal mask

No red marks on bridge of nose Under-the-nose nasal cushion, Therapy device sold separately

Highlights

Freedom of movement



By moving the tube location to the top of the head, DreamWear offers more freedom of movement during the night. Users stated they didn't need to choose a sleep position when wearing DreamWear.*

Innovative design



Prevents red marks, discomfort or irritation on the nose bridge

Unobstructed view



It's easy to wear glasses, read, watch TV, and use a computer or tablet before falling asleep.

Four-part modular design



Modular design for quick assembly and cleaning, with interchangeable cushions available in four sizes.

Soft fabric wraps



Provides a soft interface to the cheeks

Soft, flexible frame



The flexible material through which the air flows is applied gently to the skin and is pleasant to wear.

Magnetic clips

Contact your healthcare professional before you use this mask. Some medical devices may be affected by magnetic fields. The magnetic clips in this mask should be kept at least 2 in. (50mm) away from any active medical device with special attention to implanted devices such as pacemakers, defibrillators and cochlear implants. Read in-pack IFU's for more information.

HH1135/00

Specifications

Warning: Magnetic interference

- Beware with pacemakers: Read IFU's before use.

Product details

- Materials: Polycarbonate (PC) silicone, Polyester, Silicone, Spandex
- Operating pressure: 4 to 20 cmH2O
- Sound level: 18 dBA
- Number of parts: 9

Intent of use

- To be used on patients >30kg: with CPAP prescription

Cleaning: Mask

- Hand wash with mild detergent: Rinse thoroughly, air dry¹

Cleaning: Headgear, tubing

- Hand wash with mild detergent: Rinse thoroughly, hang to dry¹

Cleaning: Cushion

- Hand wash with mild detergent: Rinse thoroughly, air dry¹

Cleaning: Frame

- Hand wash with mild detergent: Rinse thoroughly, air dry¹



Issue date 2018-08-30

Version: 4.1.1

© 2018 Koninklijke Philips N.V.
All Rights reserved.

Specifications are subject to change without notice.
Trademarks are the property of Koninklijke Philips N.V.
or their respective owners.

www.philips.com

* Source: Philips User Survey 2015

** Design claim - the design of the mask is such that it does not contact the nasal bridge.

¹Please refer to Instructions For Use for more detailed cleaning instructions.